

# FEP

## MODULAR EXTERNAL FIXATOR

Lightweight  
Modular  
Sterile



# Indications

Scan the QR-Code to get more information about the device



Open, closed and intra-articular fractures.

Nonunions, limb lengthenings, deformity corrections and axial corrections.

## FEATURES

Sterile kit

Available in three sizes

Versatile

## CLINICAL ADVANTAGES

Ready for use

Wide range of application

Allows independent corrections on all planes



## Kit FEP

### Kit FEP Small - 5000060MST


Description	Qty
Body of the fixator Small	1
Self-drilling cortical screw Ø 4.0 mm L 120 mm	4
Single clamp for FEP Small and Medium	4

### Kit FEP Medium - 5000001MST

Description	Qty
Body of the fixator Medium	1
Self-drilling cortical screw Ø 5.0 mm L 140 mm	6
Single clamp for FEP Small and Medium	6

### Kit FEP Large - 5000030MST

Description	Qty
Body of the fixator Large	1
Self-drilling cortical screw Ø 6.0 mm L 180 mm	6
Single clamp for FEP Large	6

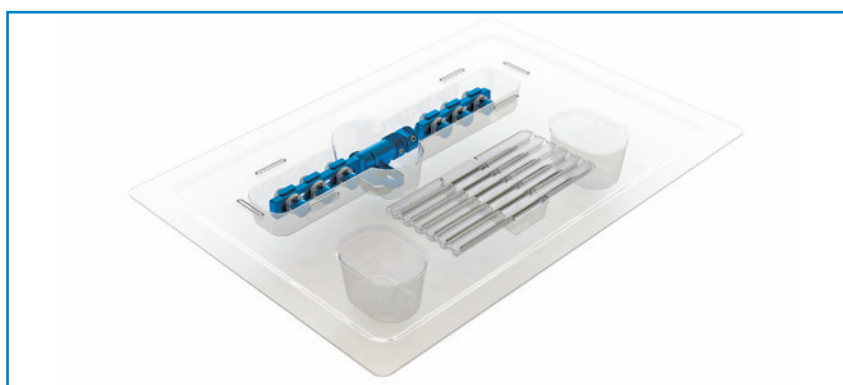

FEP kit allows independent corrections in all 3 planes. The joint unit allows a procurvatum-recurvatum correction. The compression-distraction unit allows the correction of rotational defects. The micrometric regulators allow a varus-valgus correction. Any corrections can be made independently without losing the previous correction.

The elasticity of Aluminium alloy 7012 favors the reparative processes of the fracture site while ensuring an optimal stability.



#### **HEADQUARTERS**

Via P. Gobetti, 56r - 16145 Genova  
Tel. +39 010 30801 - Fax +39 010 3080210  
[www.mikai.it](http://www.mikai.it) - [servizio.clienti@mikai.it](mailto:servizio.clienti@mikai.it)

#### **MANUFACTURING FACILITY**

Via Canestrello, 2 - 36050 Monteviale (VI)  
Tel. +39 0444 950100 - Fax +39 0444 950133  
[www.mikai.it](http://www.mikai.it) - [servizio.clienti@mikai.it](mailto:servizio.clienti@mikai.it)